

Medical Master

User manual and technical description

Model : Hospital bed table Series

Version : 2020-A-1

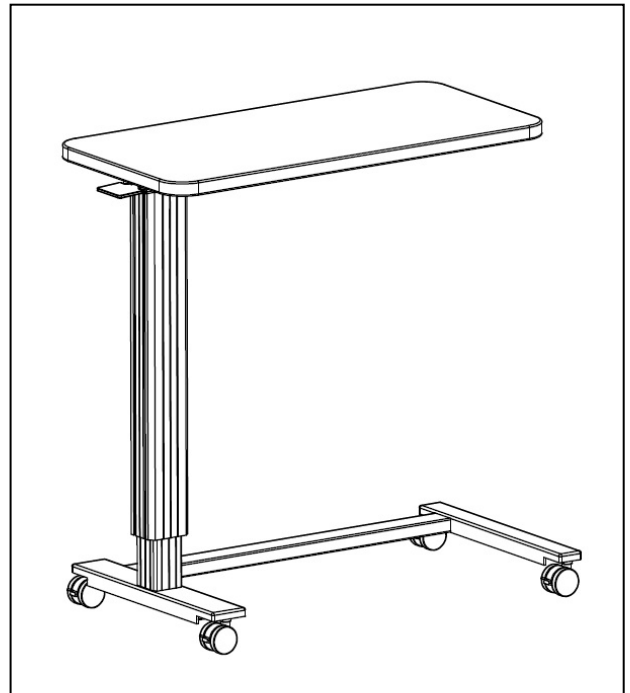


Table of Content

1 Component name and location ----- 3

 1.1 H-shaped table base ----- 3

 1.2 Π-shaped table base ----- 3

2 Hospital bed table basic specification -----

 2.1 H-shaped table base ----- 4

 2.2 Π-shaped table base ----- 5

3 Use instructions ----- 6

 3.1 Advantage ----- 6

 3.2 Operation method ----- 6

 a. Table height adjustment ----- 6


 b. Caster ----- 6

4 Consulting and service ----- 7

1 Component name and location


1.1 H-shaped table base

- 1 E1-V313 Particle Board x 1
- 2 M3.5 Phillips flat head screws x 8
- 3 Two-section lifting column with gas spring system x 1
- 4 ABS cover x 1
- 5 H shape table base x 1
- 6 Spring washer x 4
- 7 M6 Hex socket head screw x 4
- 8 2 " Caster x 4

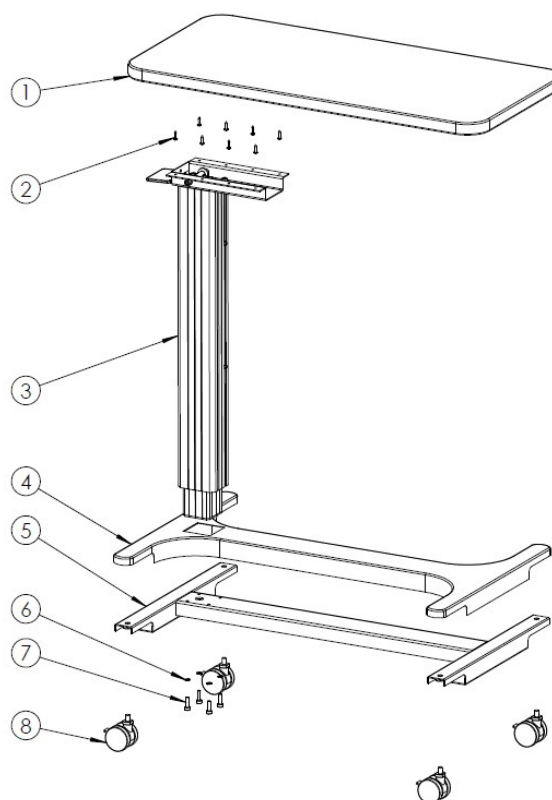


1.2 □-shaped table base

- 1 E1-V313 Particle Board x 1
- 2 M3.5 Phillips flat head screws x 8
- 3 Two-section lifting column with gas spring system x 1
- 4 □ shape table base x 1
- 5 Spring washer x 4
- 6 M6 Hex socket head screw x 4
- 7 2 " Caster x 4

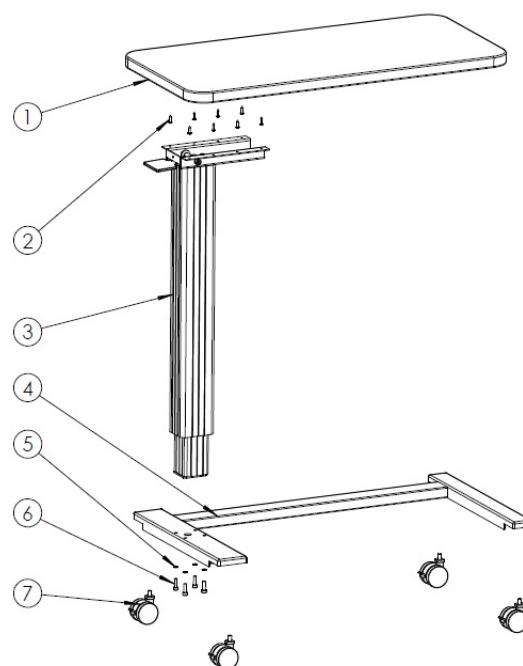


2 Hospital bed table series basic specifications

2.1 H-shaped table base

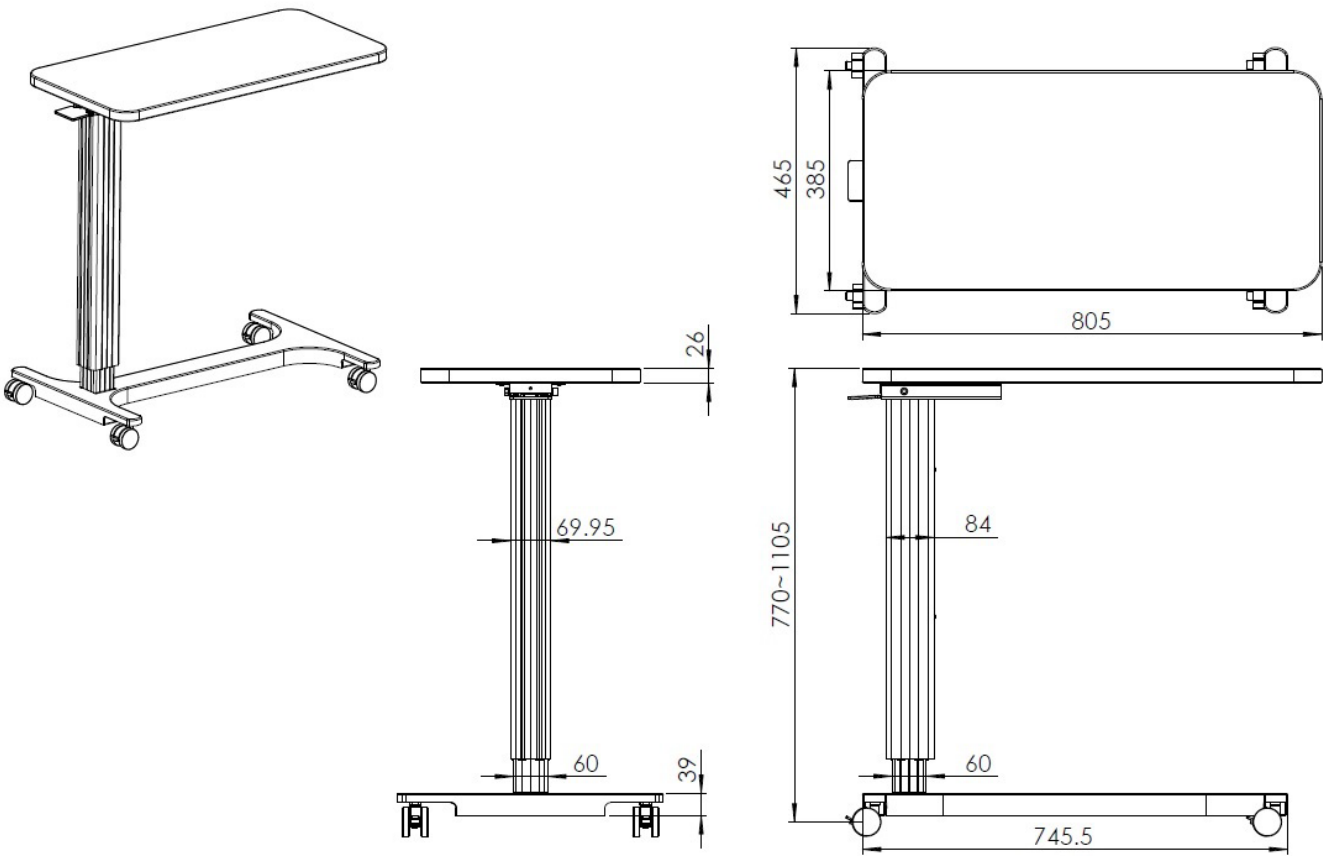Model	MDB-900
	
Specification	
Table size	L805 x W465 x H770~1105 mm
Table top	E1-V313 Particle Board
Adjustable height system	Two-section lifting column with gas spring system
Frame	Power coating / Electroplate with H shape table base
Table base cover	ABS cover
Casters	2" caster x 4
Loading capacity	Max. 20 kg

※ Tolerance : $\pm 5\%$

2.2 □-shaped table base

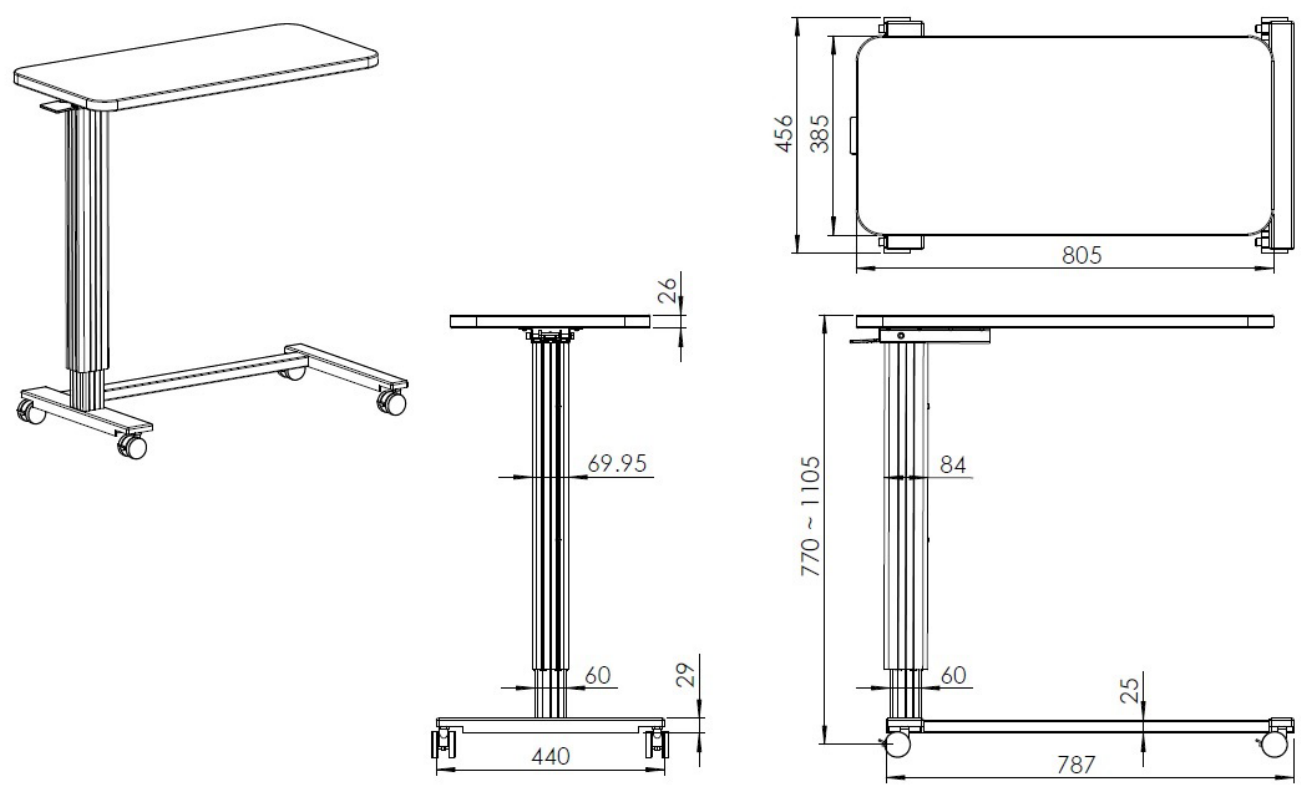Model	MDB-900	
		
Specification		
Table size	L805 x W456 x H770~1105 mm	
Table top	E1-V313 Particle Board	
Adjustable height system	Two-section lifting column with gas spring system	
Frame	Power coating / Electroplate with □ shape table base	
Casters	2" caster x 4	
Loading capacity	Max. 20 kg	

※ Tolerance : ±5%

3 Use instructions


3.1 Advantage

1. Suitable for students, mobility difficulties, single friends, use at work.
2. Pneumatic auxiliary lift, easy to operate. After the side control wrench is pulled up slightly, the table can be pressed down to adjust the height you want.
3. When caring for the elderly, place the food on the table and move it to the bedside for meals.
4. Equipped with brake wheels to facilitate users to move the table.
4. The design of □ -shaped table stand makes it easier for people with reduced mobility to use it near the table.
5. The hot-pressed table board is waterproof and drop-proof, making you happy to use.
6. Simple assembly and easy operation, can be used as a bedside reading / computer use.
7. Base ABS cover, easy to clean.

3.2 Operation method


a. Table height adjustment

Slightly pull the side control wrench up with one hand, and press down on the desktop. Feel free to adjust the height you need.

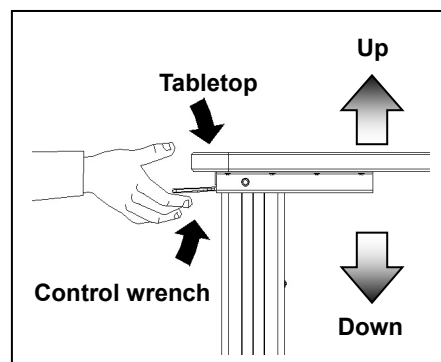

b. Caster

Press the brake pedal on the 2 "caster with your foot. When the brake pedal is down, the wheel is locked. Otherwise, the wheel is free to rotate.

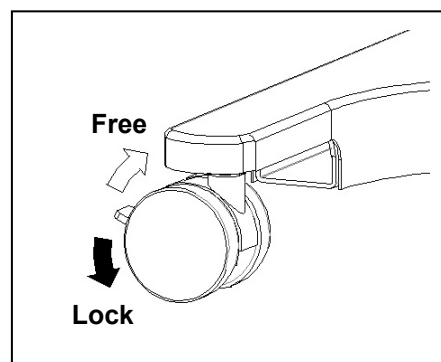

 WARNING!	<p>Risk of injury due to lifting the table too quickly!</p> <ul style="list-style-type: none"> ● Please Lift the table slowly to avoid risk of injury too quickly! ● Do not exceed the weight load of the table items ● Ensure there are no body parts or objects are between the desktop and the base.
--	---

4 Consulting and service

If you have any problems during use, please contact the following departments:

Manufacture : Medical-Msater Co., Ltd

Office

No.59-9,Ying Cheng 1st., Shi-Tung Dist, Taichung City,Taiwan (407)

Factory

No.158, Keliao Rd., Shengang Township, Changhua County,Taiwan(509)

TEL : 886-4-23171028

FAX : 886-4-23171029

E-Mail : info@medical-master.com

www.medical-master.com

Medical-Master Co., Ltd.
No. 59-5, Ying Cheng 1st, Shi-Tung Dist
Taichung City, Taiwan (40749)
Tel: +886-4-23171028