en This product contains the phthalate DEHP (applicable for tubing only). Phthalates are widely used as plasticizers in medical products and have undergone extensive testing. However, when using the device, it should be taken into consideration that possible health effects when used for long durations in pregnant women, nursing mothers and children cannot be entirely discounted. Not Made With Natural Rubber Latex.

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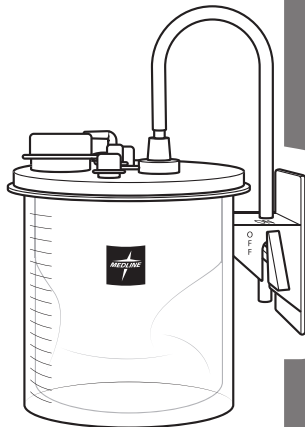
Keep Dry

Caution

Caution: Mandatory action to read the instructions for use.

Keep away from Sunlight





Medline International France SAS, 5 rue Charles Lindbergh, 44110 Châteaubriant, France

UK representative and importer in UK: Medline Industries Ltd, WA4 6HL, UK

Medline International B.V., Nieuwe Stationsstraat 10, Arnhem, 6811 KS, Netherlands


CH Rep and importer: Medline International Switzerland SARL, c/o MN & Associés SA, 1 Place de Longemalle, 1204 Genève, Switzerland

www.medline.com Medline is a registered trademark of Medline Industries, LP. Manufactured for | fabriqué pour | fabricado para: Medline Industries, LP, Three Lakes Drive, Northfield, IL 60093 USA. Made in China | Fabriqué en Chine | Fabricado en China. **1-800-MEDLINE** Distributed in Australia by: Medline International Two Australia Pty. Ltd., Sydney, Australia. V1 RH24SJP 2024-08

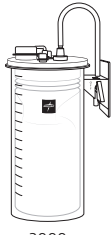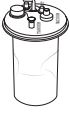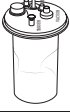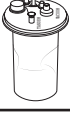
    100015003000 REF DYNDTND		
	1000 cc (mL) 1500 cc (mL) 3000 cc (mL)	REF DYNDSCL1000 (56717) DYNDSCL1500 (5711) DYNDSCL3000 (96573) OR1930
 	1000 cc (mL) 1500 cc (mL) 3000 cc (mL)	REF OR1910PG OR1920PG OR1930PG
  	1000 cc (mL) 1500 cc (mL) 3000 cc (mL)	REF OR916K OR926K OR936K
  	1000 cc (mL) 1500 cc (mL) 3000 cc (mL)	REF - OR929K OR939K
  	1000 cc (mL) 1500 cc (mL) 3000 cc (mL)	REF OR53916 OR53926 -
  	1000 cc (mL) 1500 cc (mL) 3000 cc (mL)	REF - OR53929 -
  	1000 cc (mL) 1500 cc (mL) 3000 cc (mL)	REF OR54916 OR54926 -

	1000 cc (mL) 1500 cc (mL) 3000 cc (mL)	REF DYNDCL01000 DYNDCL01500 DYNDCL03000
	1000 cc (mL) 1500 cc (mL) 3000 cc (mL)	REF DYNDCL010T DYNDCL015T DYNDCL030T

SETUP AND USE | SINGLE

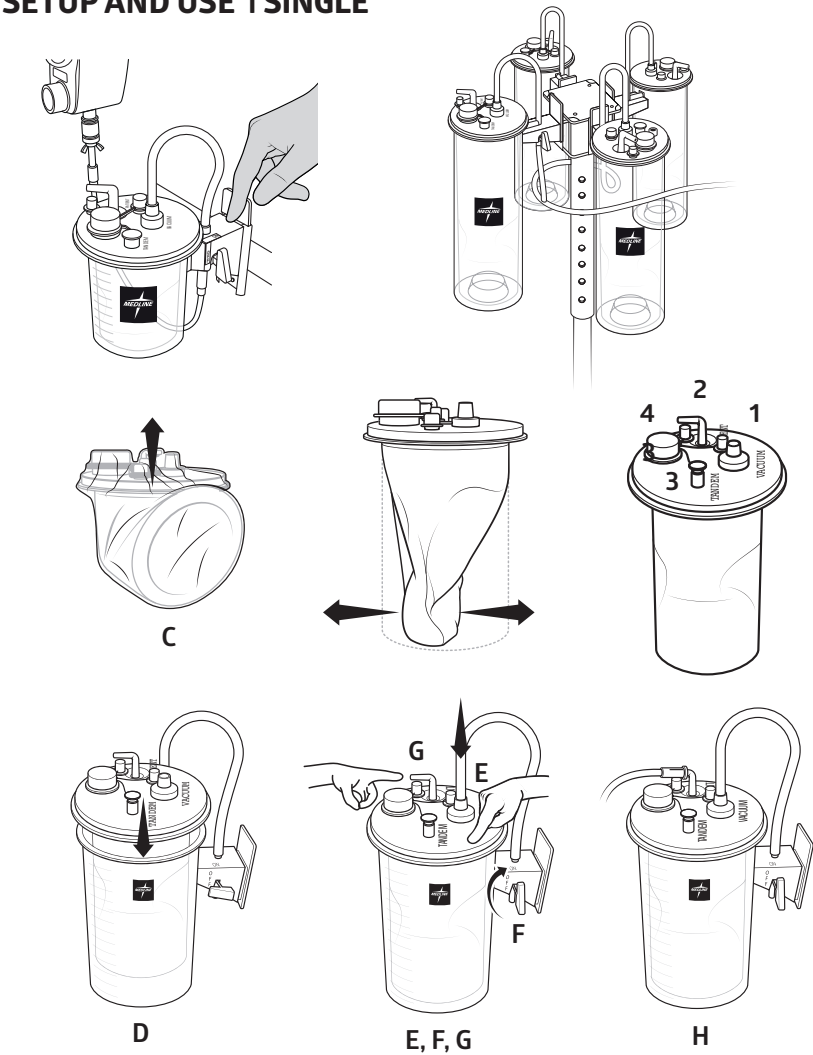SETUP AND USE | TANDEM

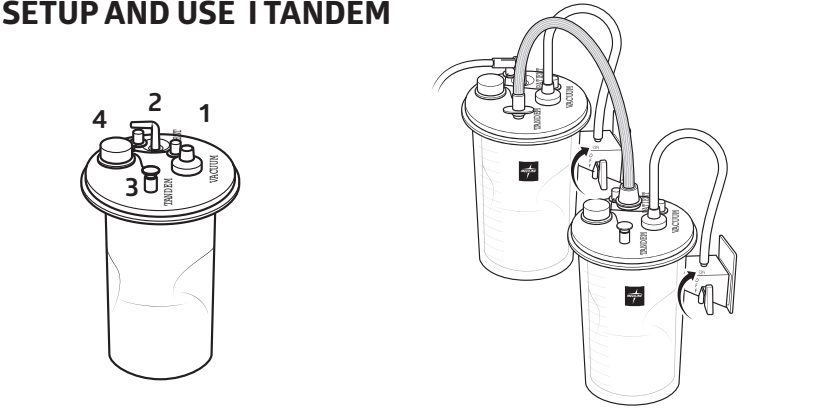en
SOFT LINERS

READ INSTRUCTIONS COMPLETELY PRIOR TO USE.

Intended Purpose:

- Medline Suction Canister System is a vacuum-powered device intended to aspirate, remove, or sample body fluids. The device is powered by an external source of vacuum.
- Intended for use on patients, within the general patient population, undergoing a procedure to sample, aspirate, or remove bodily fluids.
- Intended for use during procedures requiring the collection of fluid in a healthcare environment by a healthcare professional, who have received adequate instructions on the use of suction equipment.

Precautions:

- Avoid tilting, jarring, or sudden movements of the canister which may activate the valve and can cause premature shut-off.
- Allow for sufficient length of suction tubing (patient/liner) to prevent canister from tilting during use.
- For proper functioning of suction canisters, use appropriate vacuum tubing between canister and vacuum source that resist from kinking.
- Vacuum should only be introduced after correct device connections are made.
- Ensure vacuum suction remains on until all connecting tubing is disconnected and ports are securely capped.
- Always wear proper personal protective apparel when handling infectious waste.
- Reuse of single use devices creates a potential risk for the patient or user. It may cause contamination and/or impairment of function; which may lead to injury, illness and/or death of a patient.
- If a serious incident has occurred in relation to the device, it should be reported to Medline and the competent authority of the Member State in which the user and/or patient is established.

SETUP AND USE for Soft Liner : SINGLE

Port identification key

- 1 Vacuum Port
- 2 Patient Port
- 3 Tandem Port
- 4 Accessory Port

Directions for set up and use

- The canister should be placed in a compatible wall mount or trolley to maintain an upright position. Use of canister outside of wall mount or trolley is not recommended.
- Connect the vacuum tubing from the vacuum source to the port underneath the canister's red bracket.
- Unwrap the liner through the slit of the packaging and shake it gently to make sure the mechanic valve is not locked.
- Pull the liner down, extend it and gently insert it into the appropriately sized reusable outer canister with a pushing twist motion.
- Securely cap the tandem port and accessory port then connect the reinforced vacuum tubing from the outer canister to the vacuum port on the soft liner lid.
- Turn the vacuum 'ON' and make sure the ON/OFF valve on the outer canister is open.
- Secure the patient port through the L-shaped connector using your finger, and gently press the lid to seal the suction liner and outer canister assembly.
- Attach the patient tubing to the L-Shaped connector of the patient port.
- Make sure all unused ports are securely capped, prior to initiating suction.
- Finally, ensure that the vacuum has been created and the soft liner is fully deployed.
- The system is ready to be used.

Directions for disposal

- Once suctioning is finished, leave the vacuum 'ON' and the ON/OFF valve on outer canister open.
- Disconnect the suction tubing and L-shaped connector from the patient port. Securely cap the patient port and make sure that all other ports are securely capped as well.
- Disconnect the reinforced red vacuum tubing from the vacuum port and cap the port securely.
- Switch the valve from the outer canister into the 'OFF' position.
- Take the liner out of the outer canister.
- Dispose of the liners, according to your facility's protocols or other applicable regulations for biologically hazardous waste.

Precautions should be taken when discarding this device and disposal of the device shall be made in accordance with applicable hospital or national regulations for biologically hazardous waste.

SETUP AND USE for Soft Liners : TANDEM

Port identification

- 1 Vacuum Port
- 2 Patient Port
- 3 Tandem Port
- 4 Accessory Port

Directions for set up and use

- Make sure the canisters are in an upright position and well installed on a canister trolley.
- Connect a red manifold tubing to the port underneath each canister's red bracket, then connect the vacuum tubing from the vacuum source to the red manifold tubing underneath the holder.
- Unwrap liners through the slit of the packaging and shake it gently to unlock the mechanic valve.
- Pull the soft liners down, extend them and gently insert them into the each appropriately sized outer reusable canisters.
- Connect the red reinforced vacuum tubing from the outer canister to the vacuum port for each soft liner lid.
- Use DYNDTND tandem tubing to link all the liners as follows:
 - Starting from the first liner, insert the male connector of the tandem tube into the tandem port of the first liner and press it firmly, in order to ensure a secured seal (should not be possible to remove the male connector once installed).
 - Then, following a clockwise pattern, remove the elbow connector from the patient port of the next liner and connect it over the female connector of the tandem tube to the patient port of the second liner.
 - Repeat the tandem tube instructions (a and b) above, following a clockwise pattern until all the liners are connected.
- No tandem tube should be inserted into the tandem port of the last liner. Make sure to securely cap the tandem port of the last liner.
- Make sure that all unused ports are securely capped.
- Turn the vacuum 'ON'. Make sure the ON/OFF valve is open for all outer canisters in use.
- Secure the patient port through the L-shaped connector of the first liner, using your finger and gently press all lids to seal the suction liner and outer canister assembly.
- Attach the patient tubing to the patient port of the first liner.
- Make sure all unused ports are securely capped on each liner, prior to initiating suction.
- Finally, ensure that the vacuum has been created and the soft liners are fully deployed.
- The system is ready to be used.

Directions for disposal

- Once suctioning is finished, leave the vacuum 'ON' and the ON/OFF valve open on all outer canisters.
- Disconnect the suction tubing and the L-shaped connector from the patient port. Securely cap the patient port.
- Starting from the last liner, disconnect the female connector of the tandem tube from the patient port of the liner, and securely cap the tandem tube extremity and patient port.
- Disconnect the reinforced red vacuum tubing from each vacuum port and securely cap the ports.
- Switch valve from the outer canister into 'OFF' position.
- Take the liners out of the outer canisters.
- Dispose of the liners, according to your facility's protocols or other applicable regulations for biologically hazardous waste.

Note: Liner contents are considered potentially hazardous and proper caution should be taken when handling. For soft liners tandem setups, it is recommended that all used and unused disposable soft liners be discarded and subsequently replaced following each procedure as all contaminated by aerosols.

⚠️
WARNINGS

- Do not reuse or sterilize. If the same patient undergoes long-term treatment, the soft suction liner must be replaced at least every 24 hours.
- Do not subject the canister to active vacuum when not in use.
- Not intended for pleural/thoracic use.
- Inspect the device for the presence of cracks or damaged parts prior to use.
- Monitor vacuum gauge to assess for consistent vacuum pressure prior to and during use.
- Causes of premature shut-off include:
 - ◊ Bumping the canister or tilting at an angle can activate the shut-off mechanism inside the lid.
 - ◊ Suctioning agents that release gas, including hydrogen peroxide or povidone iodine may lead to the release of foam that can activate the shut-off mechanism or even eject fluid from the canister.

Cleaning of reusable canister

Reusable outer canister can be wiped with disinfectant.

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