



Patient Wipes

Instructions For Use

Type

Advance Patient Wipes with 2% Chlorhexidine Digluconate Antiseptic

Intended use

The wipes are designed for cleaning and reducing microorganisms on patient's body. These wipes also help in removing visible dirt and odor. Each wipe is intended for single-use only and should be used on intact skin. This product is for external use only.

This product is classified as a medical device.

Features

- Clinically assessed by dermatologists
- Perfect for pre-surgery preparations
- No need for rinsing, and free from alcohol
- Infused with Aloe Vera for skin hydration
- Versatile, resealable packaging, safe for daily body cleansing
- Thick and durable - gentle on the skin yet robust enough for all cleaning requirements
- Enriched with 2% Chlorhexidine Digluconate Antiseptic for potent antibacterial effects
- Delivers swift antimicrobial action and maintains excellent residual effectiveness for up to 6 hours
- Minimizes microbial transmission within healthcare settings

Direction For Use

1. Open the lid upwards to remove wipes.
2. Pull out one sheet and wipe desired surface for disinfection.
3. Close lid and reseal label after use to retain moisture in the pack.
4. Discard properly after use.

Warnings / Precautions

- Keep this product out of the reach of children.
- Do not use on ears, eyes, mouth and genital areas.
- Avoid contact with eyes. In case of contact with eyes, follow the provided first aid instructions.
- This product is for external use only.
- Dispose of used wipes in a waste bin; do not flush them down the toilet.

First Aid (In Case of Eye Contact):

- Hold eyes open and rinse slowly and gently with water for 15-20 minutes.
- If contact lenses are present, remove them immediately and rinse the eyes.
- Consult a doctor for treatment advice.

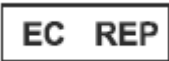
Reporting Incidents:

- Any serious incident related to the device should be reported to the manufacturer and local authorities in your country.

Storage

- Store in a cool place.
- Keep the product away from heat, direct sunlight, and open flames.

Symbol explanation

Manufacturer	Date of manufacturing	Use-by date	Batch Code	Medical device	EU representative
					
Do not reuse	Period after opening (3 month)	Keep dry	Keep away from sunlight		
					



KLEEN-PAK PRODUCTS PTE. LTD.

Address: No 2 Loyang Lane #04-03, Loyang Industrial Estate, Singapore 508913

Tel: (65) 6542 2100

E-mail: enquiries@kleen-pak.com

Brand Owner : **GLOBAL SCP MEDICAL SYSTEMS PTE. LTD.**

Address: 1 Kallang Junction, #05-01 Vanguard Campus, Singapore 339263

T : +65 6690 7330

E : idsmedglobal@idsmed.com

www.idsMED.com



Patient Wipes

Instructions For Use

Type

Comfort Patient Wipes

Intended use

The wipes are designed for cleaning and reducing microorganisms on patient's body. These wipes also help in removing visible dirt and odor. Each wipe is intended for single-use only and should be used on intact skin. This product is for external use only.

This product is classified as a medical device.

Features

- Silky, durable, and dermatologist-certied
- E-ectively antibacterial, eradicating unwelcome body odors
- Crafted without alcohol, MIT, paraben, or synthetic fragrances
- Mild and suitable for even the most delicate skin
- Easily disintegrates for hassle-free, secure disposal in most macerators
- Generously-sized for comprehensive personal hygiene, enriched with
- Aloe Vera and Vitamin E
- Enjoyable at room temperature or warmed for added comfort
- Environmentally friendly and biodegradable

Direction For Use

1. Open the lid upwards to remove wipes.
2. Pull out one sheet and wipe desired surface for disinfection.
3. Close lid and reseal label after use to retain moisture in the pack.
4. Discard properly after use.

Warnings / Precautions

- Keep this product out of the reach of children.
- Avoid contact with eyes. In case of contact with eyes, follow the provided first aid instructions.
- This product is for external use only.
- Dispose of used wipes in a waste bin; do not flush them down the toilet.

First Aid (In Case of Eye Contact):

- Hold eyes open and rinse slowly and gently with water for 15-20 minutes.
- If contact lenses are present, remove them immediately and rinse the eyes.
- Consult a doctor for treatment advice.

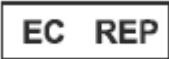
Reporting Incidents:

- Any serious incident related to the device should be reported to the manufacturer and local authorities in your country.

Storage

- Store in a cool place.
- Keep the product away from heat, direct sunlight, and open flames.

Symbol explanation

Manufacturer	Date of manufacturing	Use-by date	Batch Code	Medical device	EU representative
					
Do not reuse	Period after opening (3 month)	Keep dry	Keep away from sunlight		
					



KLEEN-PAK PRODUCTS PTE. LTD.

Address: No 2 Loyang Lane #04-03, Loyang Industrial Estate, Singapore 508913

Tel: (65) 6542 2100

E-mail: enquiries@kleen-pak.com

Brand Owner : **GLOBAL SCP MEDICAL SYSTEMS PTE. LTD.**

Address: 1 Kallang Junction, #05-01 Vanguard Campus, Singapore 339263

T : +65 6690 7330

E : idsmedglobal@idsmed.com

www.idsMED.com