HUDSON RCI®

AQUAPAK®

REF 400040,000-40F


STERILE | R

REF 400301, 400501, 400601

STERILE | A

REF 400340, 003-40
400640, 006-40

STERILE | R

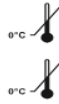Rx ONLY

Teleflex
MEDICAL

3015 Carrington Mill Blvd
Morrisville, NC 27560 USA
USA: (866) 246 6990
International: +1 919 544 8000

CE 0120 EC REP
Teleflex Medical
IDA Business and Technology Park
Dublin Rd.
Athlone, Co. Westmeath
Ireland

EN Prefilled Sterile Humidifier

Indication for use: Adds humidity to a patient's inspired breathing gases.

Instructions for use:

NOTE: Use only Hudson RCI AQUAPAK® sterile humidifier reservoirs.

1. Thread the Humidifier Adaptor onto the reservoir and tighten.
2. To open the reservoir, snap the outlet trigger upwards.

CAUTION: DO NOT TWIST.

3. Connect the Humidifier Adaptor to a flow-meter.

CAUTION: Do not over tighten wing nut.

4. Verify audible alarm function prior to use. Occlude port and listen for alarm to sound.

CAUTION: An alarm from the pressure relief valve indicates the reservoir has Accumulated too much pressure. Inspect the system for excessive blockage or a kinked tube.

5. Attach supply tubing to the reservoir outlet.
6. Turn on the flow-meter and verify the flow of gas through the device.

Warning: Single use: Do not reuse, reprocess or re-sterilize. Reuse of device creates a potential risk of serious injury and/or infection which may lead to death. Reprocessing of medical devices intended for single use only may result in degraded performance or a loss of functionality.

Contents of reservoir for inhalation therapy only, not for parenteral administration.

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