

VADI**Temperature Probe****REF VH30T150**

Length: 150cm



VADI Medical Technology Co., Ltd. Yangmei
No. 198, Ln. 298, Huandong Rd., Zhongshan Village, Yangmei Dist.,
Taoyuan City, Taiwan.
TEL: +886-3-475-9911 E-mail: sales@vadimedical.com.tw

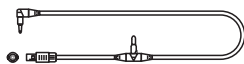
EC REP Obelis S.A.
Bd. Général Wahis 53, 1030 Brussels, Belgium.
TEL: +32-2-732-5954 E-mail: mail@obelis.net
FAX: +32-2-732-6003



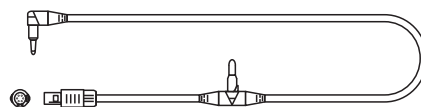
Please see the seal.

Please see the seal.

LOT Please see the seal.



Made in Taiwan

VADI**Temperature Probe****REF VH30T150****VH30T130****VH30T110**

Manual

Intended Use

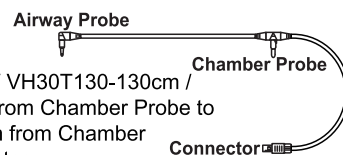
Temperature Probe is used for temperature measurements of Chamber Outlet in breathing circuit and Patient Side in Y-piece. And it transfers the measurements datum to humidifier for providing appropriate temperature and humidity.

Specification

Name: Temperature Probe

Length: VH30T150-150cm / VH30T130-130cm /
VH30T110-110cm from Chamber Probe toAirway Probe; 50cm from Chamber
Probe to the connector.

Measurement Range: 10~50°C

Tolerance: $\pm 0.5^{\circ}\text{C}$ 

Use Instruction

Please follow the following ways to install and use Temperature Probe:

1. Plug Temperature Probe into the humidifier from Connector to the sock for Temperature Probe on the humidifier.
2. Please insert Chamber Probe into the temperature line on the connector which is attached on the outlet of the humidification chamber.

※ Chamber Probe is easy-damaged. Please be careful while using Temperature Probe.

3. Please insert the airway probe into the temperature line on Y-piece.

4. Activate the humidifier (VH-3000 or equivalent device) and check the function of the temperature probe according to the manual of VH-3000 or equivalent device.

(Please refer to the manual of VH-3000 or equivalent device)

Disposition

Please DO NOT dispose Temperature Probe as normal trash. Please dispose it according to the local rule of the disposition of electronic device/accessory waste.

Clean & Sterilization

1. Please use the damp cloth applied with 70% Alcohol only.
2. Sterilizes the probe using Ethylene Oxide at 55°C before use.

Warning

- ◆ NO MODIFICATION of this device is allowed.
- ◆ Do not autoclave.
- ◆ Do not immerse other chemical disinfectants, except above recommended disinfectants in Clean & Sterilization.

VADI Medical Technology Co., Ltd. Yangmei

No. 198, Ln. 298, Huandong Rd., Zhongshan Village, Yangmei Dist., Taoyuan City, Taiwan.

TEL: +886-3-475-9911 E-mail: sales@vadimedical.com.tw

FAX: +886-3-475-9922 Website: www.vadimedical.com.tw

EC REP Obelis S.A.

Bd. Général Wahis 53, 1030 Brussels, Belgium.

TEL: +32-2-732-5954 FAX: +32-2-732-6003

E-mail: mail@obelis.net

Please see the seal.

Please see the seal.

LOT Please see the seal.

Made In Taiwan