PF-739/01/091122/02

EN

INSTRUCTIONS FOR USE

C3 EMPTYING BAG WITH SAP

REF & Description

REF	Description	Capacity	Flexi-Funnel™ Size	SAP (Fluid Absorber)	Spigot	DEHP
20460	C3 Emptying Bag w/SAP by Curas™, 2L, Dual-NRV System (incl. Spigot), Flexi-Funnel™, Cross-Hanger™, Non-DEHP, CleanPak™, WardPak™	2L	Std	✓	✓	Non-DEHP
20465	C3 Emptying Bag w/SAP by Curas™, 1.5L, Dual-NRV System (incl. Spigot), Medium Flexi-Funnel™, Cross-Hanger™, Non-DEHP, CleanPak™, WardPak™	1.5L	Medium	✓	✓	Non-DEHP
20469	C3 Emptying Bag w/SAP by Curas™, 2L, Dual-NRV System (incl. Spigot), Medium Flexi-Funnel™, Cross-Hanger™, Non-DEHP, CleanPak™, WardPak™ | PLUS	2L	Medium	✓	✓	Non-DEHP
20470	C3 Emptying Bag w/SAP by Curas™, 2L, Dual-NRV System (incl. Spigot), Flexi-Funnel™, Cross-Hanger™, DEHP Free	2L	Std	✓	✓	Non-DEHP
20461	C3 Emptying Bag, 2L, NRV, Flexi-Funnel™, Cross-Hanger™	2L	Std	X	X	Non-DEHP
20471	C3 Emptying Bag, 2L, NRV, Flexi-Funnel™, Cross-Hanger™, DEHP Free	2L	Std	X	X	Non-DEHP
20467	C3 Universal Emptying & Drainage Bag (Night Bag) w/ SAP by Curas™, 2L, 120cm, NRV, Male Connector, Cross-Hanger™, CleanPak™, Non-Sterile, WardPak™	2L	Std	✓	✓	Non-DEHP

Patented in United States of America (US11547597B2)
Additional International Patent Pending

Product Description

Disposable C3 Emptying Bag integrated with Flexi-Funnel™ for easy connection to various urine bag with bottom outlet. It comes with Cross-Hanger™ to better secure the connection. These products are latex free and are available in the following options/ features:

- 2L or 1.5L capacity
- Standard or medium size of Flexi-Funnel™
- Non-DEHP or DEHP Free material

Intended use

They are used for collecting and disposing urine from drainage bag systems. It is most suitable for collection of biohazard and Cytotoxic waste. It can also be used to collect bodily fluids. The C3 Emptying Bag helps to reduce risk of infection by providing a closed system during the urine collection. It is suitable for short-duration use and is intended for single use only.
The **SAP by Curas™** (Fluid Absorber) included to provide rapid, consistent absorption of bodily fluids and other liquids by converting these fluids to a firm, non-spill gel.

Shelf Life

5 years from date of manufacture

eIFU/UR/2046X_2047X/06102301


PF-739/01/091122/02

Caution

- Measurements are approximate volume only.
- Do not use the C3 Emptying Bag if there are any signs of tearing.
- Keep out of reach of children.
- Do not exceed the maximum liquid level.
- Do use two-hand to connect the C3 Emptying Bag with the drainable bag to ensure the connection is tight.
- Do use the Cross Hanger to secure C3 Emptying Bag onto drainable bag.

Safety Warning

- Contains or has presence of phthalate: bis 2 (2-ethylhexyl) phthalate (DEHP), except REF20470 which is DEHP Free.
- SAP warning: In case of irritation or accidental ingestion, consult physician and show product packaging. Do not eat and keep away from children.

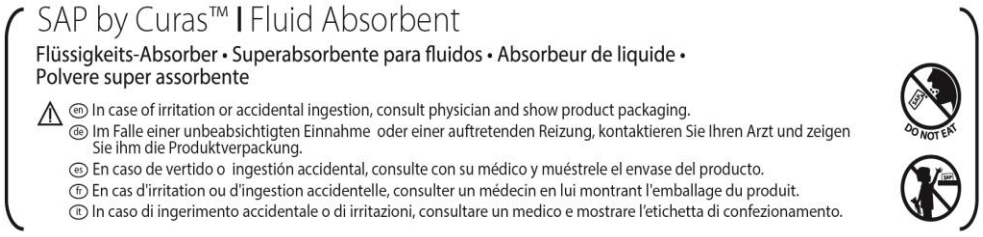

Instruction for Use (See diagram below)

1. **Preparation**
   Wash hands using soap and water, dry thoroughly and put on clean gloves before using the emptying bag.

2. **Open packaging**
   Open the product packaging (CleanPak™) and take out the emptying bag.


3. **Attaching**
   Clean the in-used urine bag's bottom outlet, and the C3 Flexi-Funnel™ using alcohol wipes. Attach the C3 Flexi-Funnel™ to the bottom outlet of drainable urine bag. Hang the Cross-Hanger™ of C3 Emptying Bag onto the bottom outlet's slider of drainable urine bag.

4. **Emptying of drainable urine bag with bottom outlet.**
   Empty the drainable bag every day or anytime when the bag is half full. Hold the bag and drain the urine from the bottom outlet into C3 Emptying Bag as follow:

   a) Drainable urine bag with Cross valve/T valve: Push the valve's slider to the other side of the housing to allow drainage. When the drainable bag is emptied, close the valve by pushing the slider back to its original position.
   b) Drainable urine bag with Lever tap valve: Push the lever down to allow drainage. When the drainable bag is emptied, return the lever to upright position.

5. **Urine Emptying Process**
   Check on the C3 Emptying Bag during urine collection to ensure the bag is not displaced when used. During urine drainage, hanging of the bag in the air with its weight supported by Cross-Hanger™ is not encouraged. Holding or placing the bag on the floor, if the length is sufficient, is encouraged. Ensure that Cross-Hanger™ is well attached to the bottom outlet of drainable bag. Check the volume content so it does not exceed the bag capacity.

6. **C3 Emptying Bag Removal**
When urine has been fully drained, close the valve of drainable bag's bottom outlet. Disconnect the C3 Emptying Bag from drainable bag then close it with Spigot. Then, clean the tip of the drainable bag's bottom outlet with alcohol wipes. Dispose the C3 Emptying Bag filled with urine immediately.

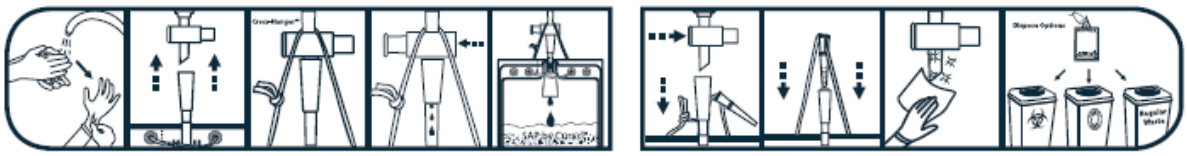

Recommended Storage Condition

- Store at dry, clean, and well-ventilated area at room temperatures between 15° to 25°C (59°- 77°F).
- Do not expose the product to sunlight, UV light, as it will change properties of the product.
- Keep away from moisture and store at dry place with relative humidity levels between 30-60% RH.

Recommended Disposal Method

- The used C3 Emptying Bag can be disposed in either a Biohazard Bin, Cytotoxic Bin or a Regular Bin, depending on the nature of the collected waste.
- Dispose used product into clinical waste bag or according to clinical waste protocol, e.g., incineration, for contaminated waste.

Product symbols

Emptying Bag Symbols



SAP by Curas™ Symbols



EC REP
Kabomed GmbH
Freiheitsstraße 43/1
A-2331 Vösendorf, Austria
office@kabomed.at

UKRP
Curas Ltd.
Country Gates, Ashton Rd.
Bristol BS3 2JH, England
info@curas.com

Curas Sdn. Bhd.
69A-B Raja Uda Business Ctr Rd 4
12300 Penang, Malaysia
info@curas.com

Curas GmbH
Oswaldgasse 6
A-1120 Vienna, Austria
info@curas.at

eIFU/UR/2046X_2047X/06102301