

flexicare
TOTAL QUALITY - TOTAL CARE

HIGH FLOW NASAL CANNULA

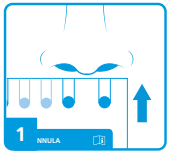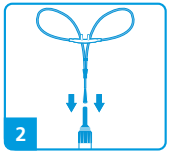Instructions for use

1. Before opening the packaging, hold the sizing guide (above) up to the nostrils to select the appropriate nasal cannula size for the patient. Caution: The cannula prongs should not completely seal the nares. If this occurs, select a smaller size.
2. Set up the breathing circuit following the manufacturer's instructions for use and attach the high flow nasal cannula to the breathing circuit limb.
3. Set the oxygen concentration and adjust the flow rate to the recommended value stated on the product labelling.
4. Ensure there is gas flow through the nasal prongs and there are no blockages/leakages in the system. Allow the gas to become warm and humidified before connecting to the patient.
5. Position the nasal prongs in the patient's nostrils and pass the tubing over the ears and under the chin. Move the sliding collar towards the chin to secure the nasal cannula.
6. Check regularly to ensure the tubing is not twisted or kinked and gas is flowing freely.

Instructions for use

- After use, dispose of in accordance with local policies and / or as advised by your local healthcare provider.
- Do not reuse this device. Reuse may result in hazards such as, but not limited to; cross contamination, damage to device and/or suboptimal performance.
- Inspect device/packaging before use. Do not use if it appears damaged.
- Any serious incident that occurs in relation to the device should be reported to the manufacturer and the competent authority of the Member State of the user and/or patient.

	PREMATURE	NEONATE	INFANT	INTERMEDIATE INFANT	PAEDIATRIC SMALL	PAEDIATRIC/ ADULT SMALL	ADULT
Max Flow (L/min)	≤8	≤8	≤8	≤8	≤20	≤40	≤40

135mm

Symbol Glossary

REF — *Catalog Number*
LOT — Batch Code
Use-by Date
MADE IN: — Country of Origin
CE 1639 — Product conforms with Directive 93/42/EEC
Do not re-use
Consult instructions for use
NON STERILE — Non-sterile
MD — *Medical device*
Do Not Use If Packaging Is Damaged
PHT DEHP — Product is not made with Phthalate DEHP
LATEX — Product is not made with natural rubber latex
Do not use blade to open
Keep Dry
Keep away from sunlight
Recycle
Manufactured by
Importer
EC REP — Authorized representative in the European Community / European Union

Flexicare Medical Limited Cynon Valley Business Park, Mountain Ash, CF45 4ER, UK
+44 (0)1443 474647 +44 (0)1443 474222 enquiries@flexicare.com


EC REP **Flexicare Deutschland GmbH** Ohlenhohnstrasse 34, 53819 Neunkirchen-Seelscheid, Germany.


IP0000-069/2
022022

125mm

flexicare
TOTAL QUALITY - TOTAL CARE
CE 1639

PRESSURE RELIEF VALVE

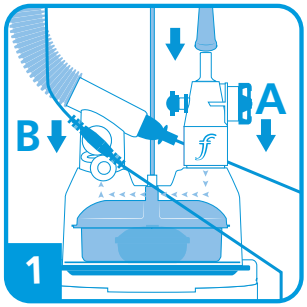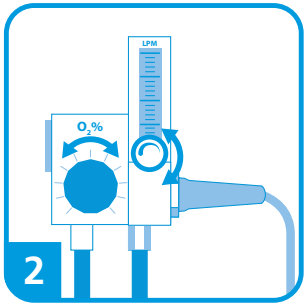Instructions for use

1. A - Position the pressure relief valve on to the humidification chamber port and connect the oxygen blender tubing.
 B - Connect the angled connector of the inspiratory limb to the other port of the humidification chamber. Set up the breathing circuit, humidification chamber and respiratory humidifier following the manufacturer's instructions.
2. Set the oxygen concentration and adjust the flow rate to the recommended value stated on the product labelling. Ensure there is gas flow through the set up and there are no blockages in the system.
3. After use, turn off the oxygen blender and humidifier, then disconnect the breathing circuit and pressure relief valve and dispose of in accordance with local policies and / or as advised by your local healthcare provider.

Do not reuse this device. Reuse may result in hazards such as, but not limited to; cross contamination, damage to device and/or suboptimal performance.


Symbol Glossary


REF Catalogue Number	LOT Batch Code	Use-by Date	Date of manufacture	CE 1639 Product conforms with Directive 93/42/EEC	Do not re-use
Consult instructions for use	NON STERILE Non-sterile	Do Not Use If Packaging Is Damaged	NO DEHP Product is not made with Phthalate DEHP	LATEX Product is not made with natural rubber latex	Do not use blade to open
Keep Dry	Keep away from sunlight	Recycle	Manufactured by		


Any serious incident that occurs in relation to the device should be reported to the manufacturer and the competent authority of the Member State of the user and/or patient.


Flexicare Medical Limited Cynon Valley Business Park, Mountain Ash, CF45 4ER, UK
+44 (0)1443 474647 +44 (0)1443 474222 enquiries@flexicare.com


IP0000-087/1 | 022020

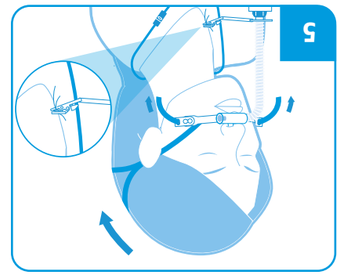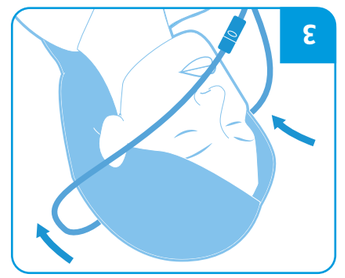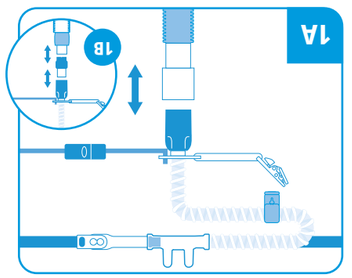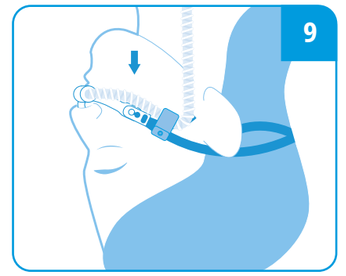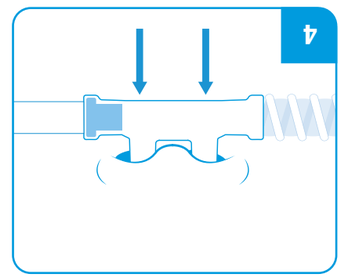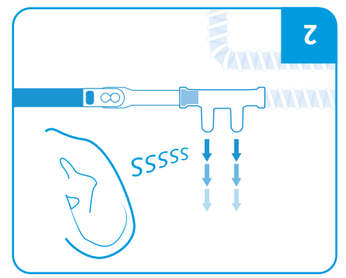
EN - HIGH FLOW NASAL CANNULA

Intended Use
Veoflo® is a high flow Nasal Cannula intended for use with a breathing system during high flow Oxygen Therapy to deliver heated and humidified gas, providing respiratory support to spontaneously breathing adult patients in the hospital.

Precautions
- Do not use if packaging is open, or if packaging or contents are damaged.
- Ensure that the Veoflo® High Flow Nasal Cannula is compatible with the breathing system. Before use, visually inspect all connections to confirm they are firmly secured and ensure there is no occlusion or leakage in the system.
- Check for gas flow from the nasal prongs before fitting the Veoflo® High Flow Nasal Cannula to the patient.
- Make sure the lanyard is fitted appropriately in order to ensure the loading of the Veoflo® High Flow Nasal Cannula is kept to a minimum.
- The lanyard is not suitable for use on patients with damaged skin.
- To prevent blockage of secretions in the patients airway, set the gas flow higher than the patient's peak inspiratory flow; this will maximise oxygen and humidity delivery.
- Check that the prongs do not completely seal the nares. If this occurs, select a smaller cannula size.
- Check regularly to ensure there are no twists in the system and gas is flowing freely.
- Check regularly for condensate in the breathing system. Drain if required.

Warnings / Cautions
- High Flow Oxygen Therapy can generate positive airway pressures. Use with care in patients where CPAP is contraindicated.
- Veoflo® High flow Nasal Cannula is single use only and intended for use up to a maximum period of 7 days.
- Always refer to manufacturers' instructions for use for all associated equipment.
- This product is single use and must not be reprocessed. This includes rinsing, washing or decontamination using gas, heat, steam, or boiling water.
- Do not use this device. Reuse may result in hazards such as, but not limited to: cross contamination, damage to device and/or suboptimal performance.
- After use, dispose of in accordance with local policies and / or as advised by your local healthcare provider.

Instructions for Use
1A. If using a Flexicare circuit, connect the cannula and connector to the heated wire breathing system; ensuring a secure fit.

1B. If required, use the Veoflo® Adapter. Slide the cannula and connector into the ridged end of the Veoflo® Adapter, then connect the smooth end of the adapter to the heated wire breathing system, ensuring a secure fit.

2. Turn on the gas supply and the heater humidifier, following the equipment manufacturers' instructions.

3. Pass the lanyard over the patient's head.

4. Place the dual elastic strap around the patient's head, positioning the nasal prongs into the nostrils. Check that the prongs do not completely seal the nares. If this occurs, select a smaller cannula size.

5. Adjust the head strap to secure Veoflo® in position. Use the clip to secure the tubing to the patient's gown to maintain the position of the system.

6. Slide the tube holder to the desired position and secure to the elastic band to avoid drag on the nasal cannula.

Any serious incident that occurs in relation to the device should be reported to the manufacturer and the competent authority of the Member State of the user and/or patient.

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148mm

Intended Use

The Veoflo® High Flow Tracheostomy Interface is intended for use with a breathing system during High Flow Therapy to deliver heated and humidified gas, providing respiratory support to spontaneously breathing adult patients in the hospital.

Precautions

- Do not use if packaging is open, or if packaging or contents are damaged.
- Ensure that the Veoflo® High Flow Tracheostomy Interface is compatible with the breathing system.
- Before use, visually inspect all connections to confirm they are firmly secured and ensure there is no occlusion or leakage in the system.
- Check for gas flow from the Veoflo® High Flow Tracheostomy Interface before attaching to the patient.
- Make sure the lanyard is fitted appropriately in order to ensure the loading of the tracheostomy tube is kept to a minimum.
- The lanyard is not suitable for use on patients with damaged skin.
- The blue deflector of the Veoflo® High Flow Tracheostomy Interface is removable.
- To prevent drying of secretions in the patient's airway, set the gas flow higher than the patient's peak inspiratory flow; this will maximise oxygen and humidity delivery.
- Check regularly to ensure there are no twists or kinks in the system and gas is flowing freely.
- Check regularly for condensate in the breathing system. Drain if required.

Part Number	Flow Range	Ambient Range	Interface Length	Tracheostomy Tube Connection
032-10-165	12–50 L/min	18–28 °C	0.3 m	15 mm Female Conical Connector
032-10-175	12–50 L/min	18–28 °C	0.3 m	15 mm Female Conical Connector

Symbols Glossary

Flexicare
TOTAL QUALITY - TOTAL CARE

Flexicare Medical Limited
Cynon Valley Business Park,
Mountain Ash, CF45 4ER UK
+44 (0)1443 474647
+44 (0)1443 474222
enquiries@flexicare.com


www.flexicare.com

IP0000-067/2 | 052020

Flexicare
TOTAL QUALITY - TOTAL CARE

veoflo®

HIGH FLOW TRACHEOSTOMY INTERFACE

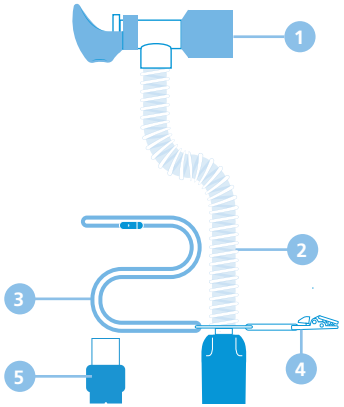1. 15mm Female Conical Connector
2. Smooth Bore Tubing **3.** Lanyard **4.** Clip **5.** Veoflo® Adaptor


Warnings / Cautions


- The Veoflo® High Flow Tracheostomy Interface is single use only and intended for use up to a maximum period of 7 days.
- Always refer to manufacturers' Instructions for Use for all associated equipment.
- Do not reuse this device. Reuse may result in hazards such as, but not limited to; cross contamination, damage to device and/or suboptimal performance.
- This product is single use and must not be reprocessed. This includes rinsing, washing or decontamination using gas, heat, steam, or boiling water.
- Do not obstruct exhalation ports.
- Do not use with an oxygen / air entrainer.
- After use, dispose of in accordance with local policies and / or as advised by your local healthcare provider.

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Instructions for Use

1. Set up the breathing circuit, humidification chamber and respiratory heater humidifier following the manufacturers' instructions.
 - Place the lanyard around the patient's neck or tape to the shoulder to reduce tube drag.
2. Connect the bottom end of the tracheostomy interface to the breathing circuit, ensuring a secure fit.
3. Turn on the gas supply and respiratory heater humidifier, following the manufacturers' instructions.
4. Confirm gas flow and allow the gas to become warm and humidified before attaching the patient end of the tracheostomy interface to the tracheostomy tube.
5. Use the clip to secure the tubing to the patient's gown to support the weight of the system.

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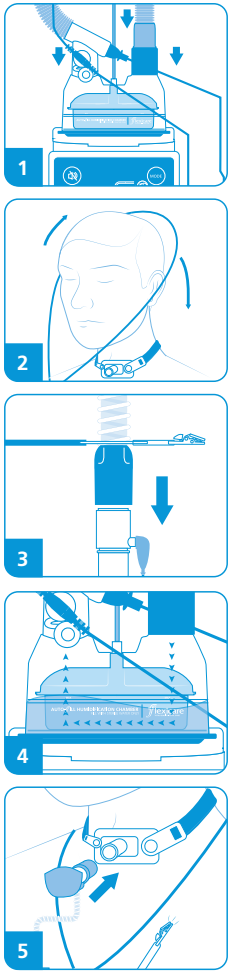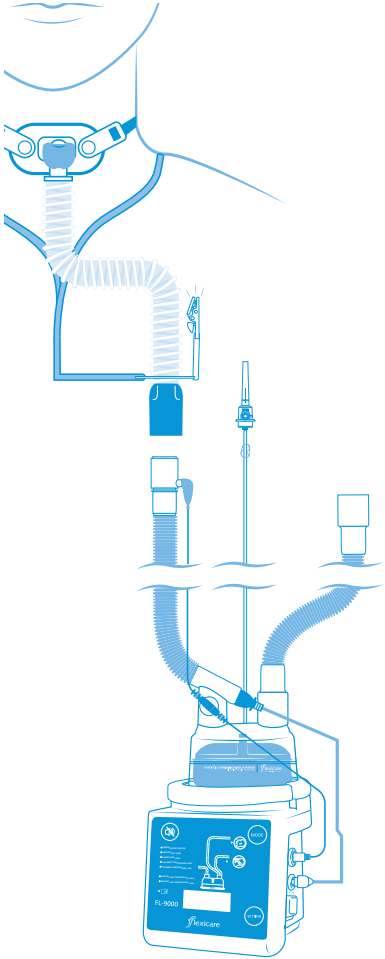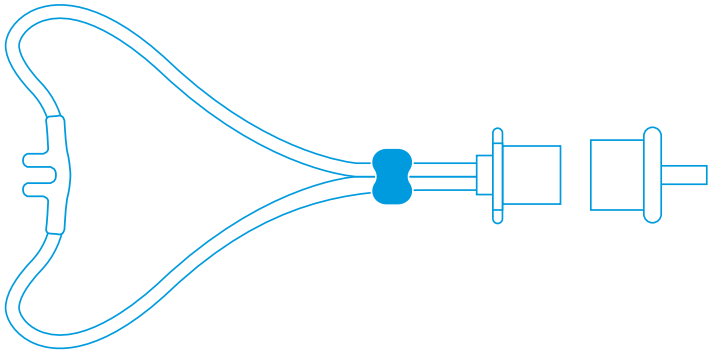
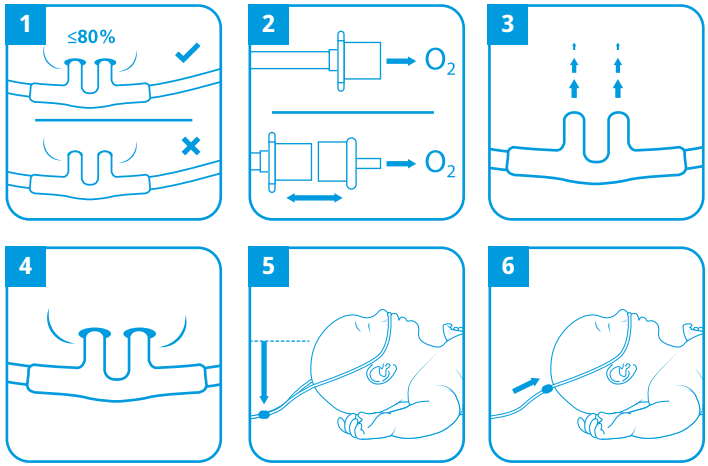
INTENDED USE

NioFlo nasal cannula is intended for use for spontaneously breathing neonate/infant and paediatric patients in the hospital who require supplementary oxygen or respiratory support.

It is used to provide oxygen with the following methods:
1. Heated or unheated, humidified low flow oxygen therapy
2. Heated, humidified high flow oxygen therapy

PRECAUTIONS

- Do not use if packaging is open, or if packaging or contents are damaged.
- Ensure the correct setup is chosen. An optional oxygen tubing adapter is included.
- Before use, visually inspect all connections to confirm they are firmly secured and ensure there is no occlusion or leakage in the setup.
- Check for gas flow from the nasal prongs before fitting the NioFlo nasal cannula to the patient.
- Check that the prongs do not completely seal the nares. If this occurs, select a smaller cannula size.
- Check regularly to ensure there are no twists or kinks in the system and gas is flowing freely.
- Replace NioFlo nasal cannula every 14 days, or less if dictated by local policies and guidelines.

WARNINGS / CAUTIONS

- Always refer to manufacturers' Instructions for Use for all associated equipment.
- This product is single use and must not be reprocessed. This includes rinsing, washing or decontamination using gas, heat, steam, or boiling water.
- Do not reuse this device. Reuse may result in hazards such as, but not limited to; cross contamination, damage to device and/or sub-optimal performance.
- After use, this product is contaminated waste and may be a potential biohazard. Handle and dispose of in accordance with hospital policy and applicable local guidelines and regulations.
- Do not position the tubing underneath the patient's head.

MICRO PREMATURE
OD 2.2mm

PREMATURE
OD 2.2mm

NEONATAL
OD 3.5mm

INFANT
OD 4.0 mm

SMALL
OD 5.3 mm

MEDIUM
OD 5.7 mm

LARGE
OD 6.5 mm

Symbol Glossary

Catalogue number

Batch code

Use by date

Product conforms with Directive 93/42/EEC

Do not re-use

Consult instructions for use

Non-sterile

Do not use if packaging is damaged

Product is not made with Phthalate DEHP

Product is not made with natural rubber latex

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Flexicare Medical Limited
Cynon Valley Business Park,
Mountain Ash, CF45 4ER UK
+44 (0)1443 474647
+44 (0)1443 474222
enquiries@flexicare.com

www.flexicare.com

IP0000-058/1 | 032020